


BIOSTAR[®] 4-Aligner

Increase your productivity without compromising on quality.



A new level in pressure moulding technology.

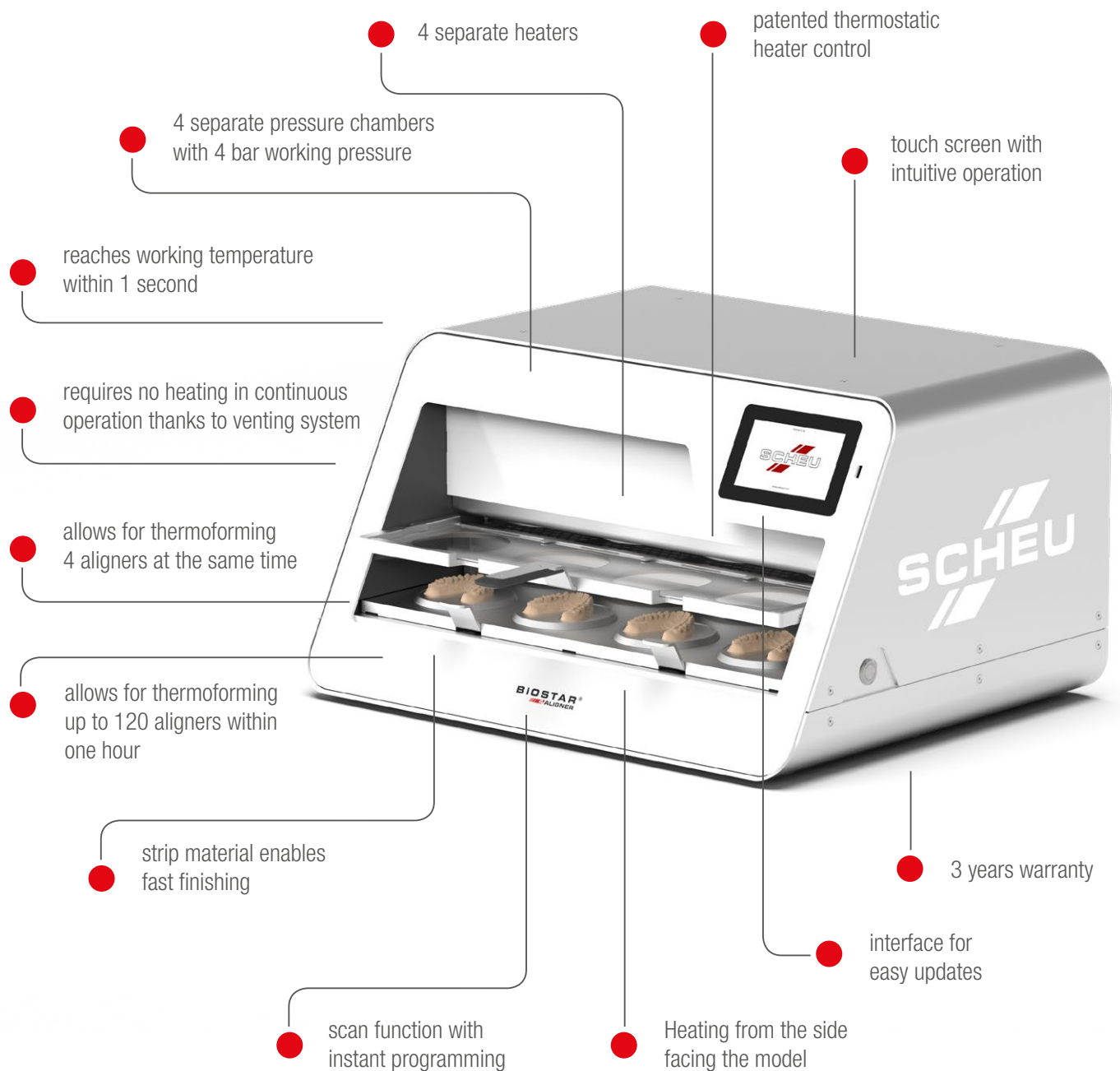

A win for your daily laboratory and practice routine!

Productive, precise and flexible - start your own high-end production now, enabling you to thermoform up to 4 aligners at the same time.



High productivity

Make your workflow more comfortable and efficient. With the BIOSTAR® 4-Aligner it is possible to thermoform up to 4 foils at the same time and thus 120 aligners in one hour.



Precise results

Each foil is individually heated in its own pressure chamber and thermoformed at 4 bar. The foils fit perfectly to the tooth rim, thus achieving optimum pressure moulding results with the usual high aligner quality.



Flexible processes

With the BIOSTAR® 4-Aligner, you can choose from various SCHEU-DENTAL foil formats depending on your individual preference: Round foils (Ø 120-125 mm) or rectangular stripes exclusively available from SCHEU-DENTAL for quadruple application (499 x 114 mm).

Features BIOSTAR® 4-Aligner

| | |
|---------------------------|--------------------|
| BIOSTAR® 4-Aligner | 95-230 V |
| REF | #3030 |
| Power | 1450 watt |
| Working pressure | 4 bar |
| Dimensions (W x H x D) | 612 x 337 x 524 mm |
| Weight | 55kg |

An unbeatable combination.

Perfect results with the innovative foil format from SCHEU-DENTAL.



The stripes, which are exclusively available from SCHEU-DENTAL, are designed for quadruple application and are provided with a scannable code. As with all SCHEU-DENTAL thermoforming units, all you have to do is scan the code on the stripes and the parameters for the thermoforming process, such as heating and cooling times, are preset. Round foils (Ø 120-125 mm) can also be thermoformed easily with the BIOSTAR® 4-Aligner to produce up to 2 aligners at the same time.

Best results for different indications:



Aligner

Material:

DURAN®+, CA® Pro+*, UNIQ pd, UNIQ+ pd

* not available in the USA

Of course, our service does not end when you purchase one of our devices.

If you have any questions or experience any problems in your daily working routine, please feel free to contact the SCHEU-DENTAL team. You can call us or arrange an appointment with our field sales staff to discuss specific questions.

